




RE-USABLE ISOLATION GOWN

(Maximum 75 wash cycles)

PRODUCT CODE:

Size Large: **222-L**

Size Extra Large: **222-XL**

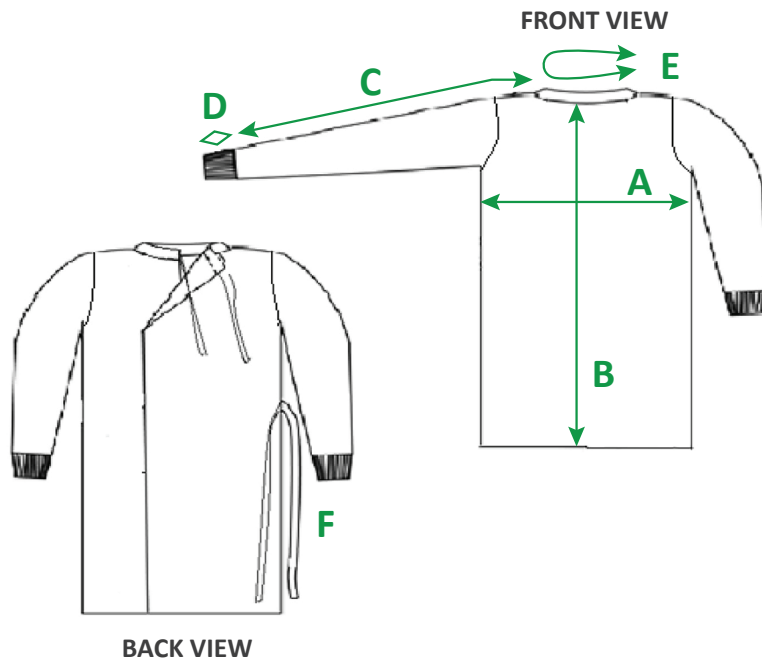

MEASUREMENTS ARE IN CM

Acceptable tolerance = 2cm

GRADE RULE

CM	MEASURING POSITION	L	XL
A	1/2 chest at under arm	70	79
B	Gown length at Centre Back including binding	119	123
C	Sleeve length - top seam from shoulder neck point to cuff (excluding cuff)	78	80
D	Cuff depth	6	6
E	Binding measurement at neck line 37mm	115	120
F	Belt Measurements 37mm	200	220
Neck binding colour		GREEN	BLUE


COLOUR OPTIONS




Yellow



Dark blue



Light blue




Lilac

Note: Colour block are an indication only and fabric actual will vary

BROGO01

WASHING AND DRYING INSTRUCTIONS

MICROFIBRE FABRIC ISOLATION GOWN (ISG)

WASHING PROCEDURE	
Water temp	38 – 65 °C
Cycle Time	36 min – 45 min
Chemicals Used	Low foaming liquid laundry detergent Alkaline Booster Oxygen Bleach Sour/Neutralizer, but not Softener
Tumble Dry	Medium heat (55 - 70°C)
	DO NOT IRON

ADDITIONAL NOTES:

- Isolation gowns are validated to 75 wash cycles;
- Above mentioned wash instructions adhere to the SANS10146 Laundry standards. This states that all hospital linen should be disinfected during the cleaning process; either a thermal or chemical disinfection. The lower the wash temperature the longer the contact time.
- The Isolation gowns will “crinkle” if washed and dried at too high temperatures